

DEFIGARD TOUCH 7 PHYSIOGARD TOUCH 7

Defibrillator / Monitor

Maximum performance, minimum weight, touch screen control: the power to save lives



SCHILLER
The Art of Saving Lives

DEFIGARD TOUCH 7 PHYSIOGARD TOUCH 7

A compact combination of state-of-the-art defibrillation technology and comprehensive monitoring functions

The DEFIGARD Touch 7 and PHYSIOGARD Touch 7 are indispensable tools for rescuers:

- available as defibrillator/monitor or only monitor
- extremely lightweight and compact
- high-end touch screen
- comprehensive monitoring functions
- equipped with enhanced data transmission technology



**DEFIGARD TOUCH 7
MONITOR / DEFIBRILLATOR**



**PHYSIOGARD TOUCH 7
MONITOR**

UNPARALLELED TOUCHSCREEN

The DEFIGARD Touch 7/PHYSIOGARD Touch 7 (DGT7/PGT7) is the first emergency monitor/defibrillator with a touchscreen. Thanks to Projected Capacitive Touch technology, the display enables highly responsive operation, even when wearing gloves. The software environment Touch'n'Save offers unmatched user-friendliness. The 7-inch touchscreen is particularly bright. Its high resolution and the exceptional display quality enable optimal display readability at all times.

CAN BE USED ANYWHERE

Thanks to SCHILLER's extensive experience in the rescue sector, the DGT7/PGT7 has been developed so that it can also be used in extreme conditions. It is the ideal companion for all land and water interventions.

Through the lightweight structure, it also offers new possibilities for air rescue operations. Its low weight ensures significant fuel saving when used on an air plane or a helicopter.

With its mounting bracket – which also serves as battery charger – the DGT7/PGT7 is resistant to shocks and vibrations, thus complying with land and air rescue operation standards.

CONNECTED AT ALL TIMES

The device can perform high speed data transmission at all times. Depending on the location, transmission can be performed using Wi-Fi or cellular network.

The safety standards WPA2 Enterprise guarantee secure data protection.

PRACTICAL AND SOPHISTICATED

In the functional carrier bag, each accessory has its place, is easily accessible and can readily be used, as each part is already connected. Moreover, the device can be tilted for optimal convenience.


OPERATIONAL AT ALL TIMES

The DGT7/PGT7 is ready to use at all times, thanks to regularly performed self-tests and to a RFID tag that ensures the monitoring of the defibrillator electrodes' expiry date.



**ONLY 2.4 KG.
LIGHTWEIGHT, COMPACT
AND USER-FRIENDLY:
YOUR PERFECT COMPANION!**

ADAPTIVE USER INTERFACE

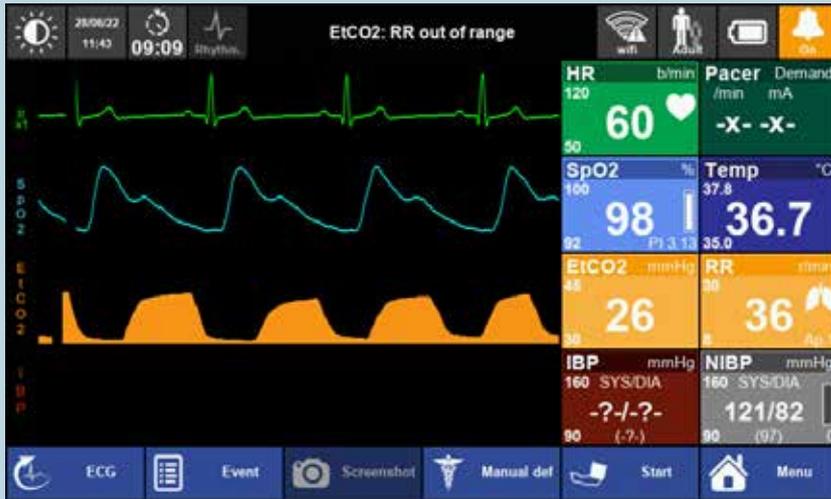

PHYSICIAN MODE This display, showing values and waveforms as well as a manual defibrillation function, gives the physician a comprehensive overview of the patient.



RESCUE MODE With a simple swipe gesture, the user can move from the physician display to the rescue display. The values are displayed numerically and are well visible, even from a distance. Moreover, a semi-automatic defibrillation function is also available for all users.

ECG

The DGT7 / PGT7 enables continuous ECG monitoring as well as 12-lead ECG diagnostic recording and ETM interpretation.

CAPNOGRAPHY

The EtCO₂ value is a parameter to check blood circulation and deduct the efficiency of CPR.

The DGT7 / PGT7 offers Mainstream or Sidestream measurement using Masimo sensors to display EtCO₂, Respiration Rate and Apnea information.

SPO₂, SPCO AND SPMET

The Masimo Rainbow technology for SpO₂, SpCO and SpMet measurements is available as an option.

The carboxyhaemoglobin (SpCO) enables the non-invasive and immediate detection of heightened levels of carbon monoxide in the blood. This ensures rapid diagnosis and treatment of patients with carbon monoxide poisoning.

The SpMet measurement enables the non-invasive and immediate detection of high levels of methemoglobin in the blood.

NON-INVASIVE BLOOD PRESSURE MEASUREMENT

Arterial blood pressure measurement in children and adults, as well as a venous block function.

INVASIVE BLOOD PRESSURE MEASUREMENT

Arterial or venous invasive blood pressure measurement with various pressure sensors.

TEMPERATURE

Continuous measurement of the temperature through skin sensors or endocavity probes.

DIRECT DATA TRANSMISSION



INTERVENTION SITE

Recording of a 12-channel ECG on site and immediate transmission to the referring physician.




MEDICAL REGULATION

Reading and analysis of the data transmitted from the intervention site.



EMERGENCY TRANSPORT

The emergency doctor can implement the instructions immediately and thus help saving essential minutes.

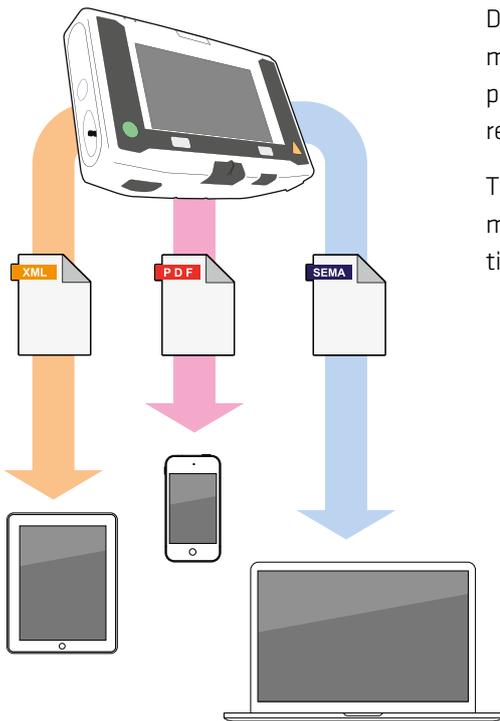

Direct data transmission ensures optimal care of the patient. The rescuer or physician can remain in contact with the regulating doctor.

The DGT7 / PGT7 enables data transmission during and after the intervention.

DURING THE INTERVENTION

- ❖ 12-lead ECG with ETM interpretation
- ❖ Trends
- ❖ Screenshots
- ❖ PDF reports
- ❖ Connection with ePCR

AFTER THE INTERVENTION

- ❖ Comprehensive intervention data including 12-lead ECG with ETM interpretation and SpO₂ curves, events and trends.



CONNECTED AT ALL TIMES

- ❖ This device is compatible with the SCHILLER Update Server, which makes the updating through the network easier than ever before.
- ❖ Manage the entire configuration of the DGT7 / PGT7: the devices can be organised into groups, build configuration templates and be deployed easily to many devices at once.
- ❖ Device self-test results are automatically sent by e-mail, for a better pool management and immediate action in case of error.

- ❖ Thanks to its Bluetooth interface, the DGT7 / PGT7 offers the possibility to print Resting ECGs, live ECGs sequences, etc. wirelessly on a thermal printer.



DEFIGARD TOUCH 7



ARGUS LIFEPOINT: The reusable CPR-Feedback-Sensor is easy to position.



ASSISTANCE IN PROVIDING CARDIAC MASSAGE

Thanks to the ARGUS LifePoint CPR feedback sensor, cardiac massage monitoring is made easy. The following elements are displayed with various colors. Available information:

- Compression depth and rate
- Pressure relief at the end of each compression¹



MANUAL OR SEMI-AUTOMATIC DEFIBRILLATOR WITH UP TO 200 JOULES

The DEFIGARD Touch 7 can be operated in semi-automatic or manual mode (synchronously or asynchronously) with Multipulse Biowave (biphasic), which is achieved with exceptional efficacy and safety, as well as with extremely low energy input.

This makes this device ideal for all patient types. It identifies paediatric electrodes and adjusts the emitted energy automatically.



EXTERNAL PACEMAKER

- Transcutaneous stimulation via pads in Fix and Demand mode
- Wide range of settings

¹ Not valid for Great Britain, Germany and France

PROVIDES ALL FUNCTIONS²



AED
MANUAL DEFIBRILLATOR
CPR FEEDBACK
PACEMAKER
12 LEADS ECG
SPO₂, SPCO, SPMET
NIBP
ETCO₂ MAINSTREAM
ETCO₂ SIDESTREAM
IBP
TEMPERATURE
WIFI
CELLULAR NETWORK
BLUETOOTH
USB
COMPACT
LIGHTWEIGHT
TOUCH SCREEN
HIGH CONTRAST DISPLAY
NUMERIC VALUE VIEW
AIRCRAFT COMPATIBLE
AMBULANCE COMPATIBLE

² To configure your perfect rescue unit, please contact your SCHILLER sales representative.



SYNOVO MED

MEDICAL SOLUTIONS

Exclusive distributor of Schiller solutions in Canada:
discover cutting-edge, innovative, and
reliable medical technologies.

Our team is here to help you find
the solutions best suited to your needs.

CONTACT US!

SYNOVO GROUP NORTH AMERICA

schiller-ca@synovo-group.com

(819) 539-1825

www.synovomed.com

1250 Av. de la Station, Shawinigan, QC G9N 8K9, Canada